

## EXTERNAL PROBE ADAPTER

CAP-MP

EXTERNAL ADAPTER FOR  
CONTINUOUS VALIDATION PROCESSES

ACCESSORY



# External probe adapter


## APPLICATION

- + It allows access for an external probe to obtain a temperature reading from the culture medium, independent of the reading taken by the equipment's flexible probe.




## CHARACTERISTICS

- The port is located on the door of the media preparator.
- The TP-VAL-MP temperature probe, compatible for this port, is available upon request.
- The external probe adapter is composed of three parts:

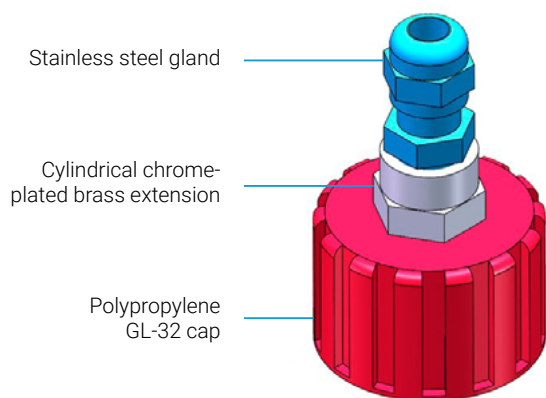


## INSTALLATION

1. Unscrew one of the red caps.\*

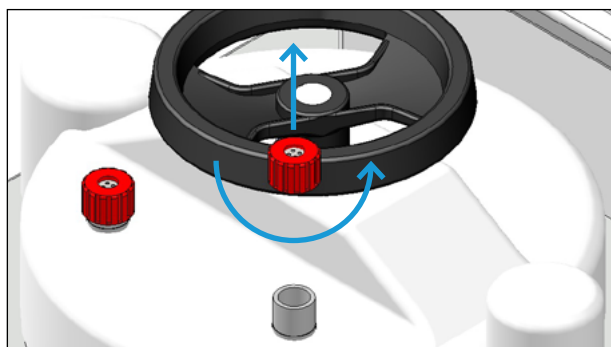

2. Tighten the CAP-MP accessory well.

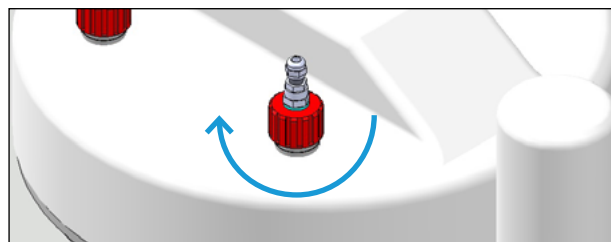

## Accessory



**EXTERNAL PROBE FOR QUALIFICATION**  
Ref. TP-VAL-MP

## Technical data

Reference	CAP-MP
Materials	Stainless steel, chrome-plated brass, and polypropylene
Height mm	65
Diameter mm	40
Compatible probe diameter mm	3-5

\*The external probe adapter is compatible with both access points; however, we recommend installation in the access point located on the right, as shown in the images.

