

## MAGNETIC STIRRING PADDLE SYSTEM

### ULTRA-STIRR

SPECIAL ADAPTATION OF THE  
MAGNETIC STIRRER WITH EXPANDED  
PADDLES

ACCESSORY



# Magnetic stirring paddle system

## APPLICATION

- + Tangential flow paddle agitation system, designed to prevent the generation of turbulence, is ideal for achieving uniform and efficient thermal exchange for the entire solution, as well as providing a gentle treatment of the product.




Must be installed at our factory.

## RECOMMENDED FOR:

- + Applications that require processing high viscosity solutions such as aqueous solutions with starch, oats, or other vegetable flours.

## MAINTENANCE

- Periodic maintenance is efficiently integrated with the cleaning cycle itself, ensuring that all components of the media preparator maintain their optimal performance.
- Additionally, to ensure uninterrupted performance and the good condition of the system, it is recommended to remove the paddle system once a week and carry out a deeper cleaning.

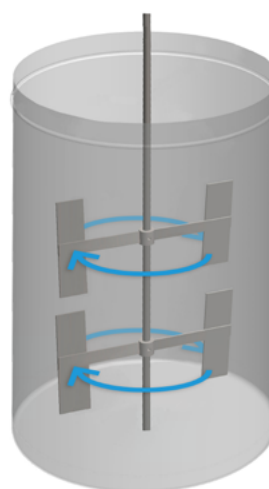

## Technical data

Reference	ULTRA-STIRR-40	ULTRA-STIRR-60	ULTRA-STIRR-80	ULTRA-STIRR-100
Dimensions Ø x H mm	270 x 500	270 x 730	300 x 640	300 x 780
Compatible media preparator	AE-40-MP	AE-60-MP	AE-80-MP	AE-100-MP

